

# Product Guide - Electrodes



**AXELGAARD**  
MANUFACTURING CO., LTD.

# Our Mission:

## Innovation and Quality

Axelgaard Manufacturing's mission is to deliver innovative and high quality electrodes to clinicians and patients worldwide.



Since 1985, Axelgaard has proven its reputation as a global leader in the design, development, and manufacturing of electrodes. Axelgaard products successfully improve a wide range of conditions for patients and athletes.

Our advanced design process is patient focused, utilizing patented technology to promote comfort, safety, and treatment outcome.

Axelgaard provides value to customers through precise and state-of-the-art manufacturing technology, consistently resulting in the highest quality products.



We invite you to learn more, connect and collaborate:  
**[www.axelgaard.com](http://www.axelgaard.com) | [support@axelgaard.com](mailto:support@axelgaard.com)**

Our products are patent protected, see:  
**[www.axelgaard.com/PatentList](http://www.axelgaard.com/PatentList)**

# Index

## Neurostimulation Electrodes

|                                       |    |
|---------------------------------------|----|
| Educational Videos .....              | 4  |
| MultiStick® Hydrogel Technology ..... | 5  |
| PALS® .....                           | 6  |
| PALS® Foam .....                      | 7  |
| PALS® StimTrove® .....                | 7  |
| PALS® Blue .....                      | 7  |
| UltraStim® .....                      | 8  |
| UltraStim® Wire .....                 | 10 |
| UltraStim® Dual .....                 | 10 |
| UltraStim® X .....                    | 11 |
| UltraStim® X Blue .....               | 11 |
| UltraStim® Snap .....                 | 12 |
| UltraStim® Garment .....              | 13 |
| ValuTrove® Cloth .....                | 14 |
| ValuTrove® Foam .....                 | 14 |
| ValuTrove® X .....                    | 15 |
| ValuTrove® Lite .....                 | 15 |
| EconoStick® Hydrogel Technology ....  | 16 |
| EconoStim® .....                      | 17 |

## Iontophoresis Electrodes





|                 |    |
|-----------------|----|
| Ionto480® ..... | 18 |
|-----------------|----|

## ECG and EMG Electrodes

|                    |    |
|--------------------|----|
| Little PALS® ..... | 19 |
|--------------------|----|

# Educational Videos

Electrode Placement & Functional Movement

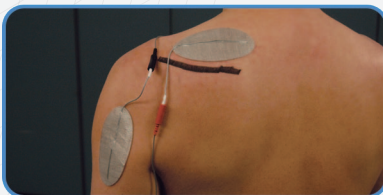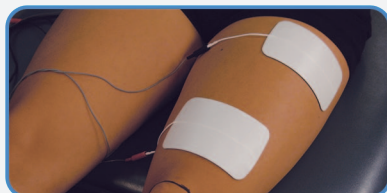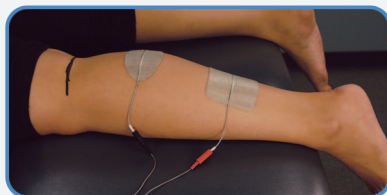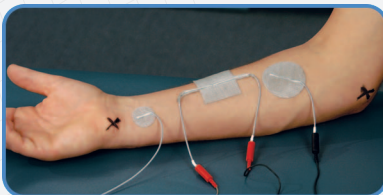

We are pleased to sponsor this educational video series presented by Dr. Lucinda Baker, Associate Professor at USC Division of Biokinesiology and Physical Therapy. These videos provide comprehensive information on the preparation and use of electrodes for effective neurostimulation treatment.

All electrodes used in the video presentation are our patented PALS® neurostimulation electrodes.

## Multiple Viewing Options


We trust these videos will assist efforts in providing the most effective treatment to your patients. DVD orders and **free** viewing/download are available on our website.

View Online at [www.axelgaard.com/education](http://www.axelgaard.com/education)

# MultiStick<sup>®</sup>

## Patented Multi-layer Adhesive Hydrogel

Axelgaard electrodes feature our multi-layered patented MultiStick hydrogel. The formulation of the hydrogel top layer is optimized for bonding to the electrode to prevent delamination. Controlled scrim placement provides consistent structural support. And the bottom layer is specifically formulated for multiple applications to skin.

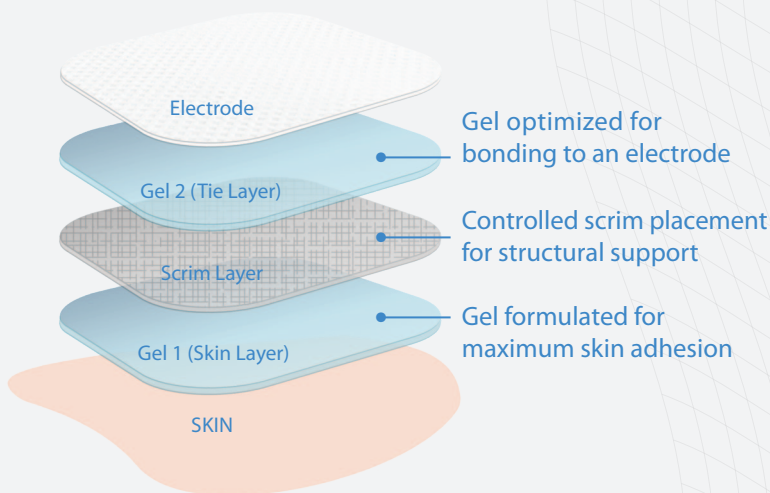


### Superior Gel Adhesion and Performance

This multi-layer approach ensures the highest quality electrodes in the industry and avoids the delamination and poor adhesion common in single layer gels.



# PALS<sup>®</sup> Neurostimulation Electrodes with MultiStick<sup>®</sup> Hydrogel

Designed for all FES and TENS applications

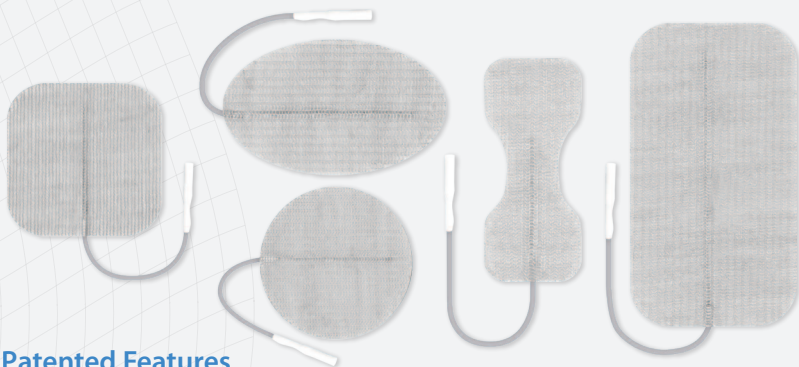

## Patented Features

- Patented stainless steel conductive fabric provides the ultimate in flexibility and conformity.
- Patented dual-layer MultiStick hydrogel provides comfort and optimizes multiple applications to the skin.

## Optimum Current Distribution

- Enhances patient comfort during stimulation.
- Disperses current evenly across the electrode.
- Eliminates stinging, edge-biting, and hot spots.


| PART # | SHAPE                                                                               | DESCRIPTION                 | #ELC/PK | #PK/CS | #ELC/CS |
|--------|-------------------------------------------------------------------------------------|-----------------------------|---------|--------|---------|
| J10R00 |    | 1.0" (2.5cm)                | 4       | 10     | 40      |
| 879100 |   | 1.25" (3.2cm)               | 4       | 10     | 40      |
| 879200 |  | 2.0" (5cm)                  | 4       | 10     | 40      |
| 879300 |  | 2.75" (7cm)                 | 4       | 10     | 40      |
| 891100 |  | 1.1" x 2.5" (2.9cm x 6.4cm) | 4       | 10     | 40      |
| 891200 |  | 1.3" x 2.1" (3.3cm x 5.3cm) | 4       | 10     | 40      |
| 895220 |  | 2.0" x 2.0" (5cm x 5cm)     | 4       | 10     | 40      |
| 895240 |  | 2.0" x 3.5" (5cm x 9cm)     | 4       | 10     | 40      |
| 895402 |  | 2.0" x 4.0" (5cm x 10cm)    | 2       | 10     | 20      |
| 895250 |  | 2.0" x 5.0" (5cm x 13cm)    | 2       | 10     | 20      |
| 895340 |  | 3.0" x 4.0" (7.5cm x 10cm)  | 2       | 10     | 20      |
| 896230 |  | 1.5" x 2.5" (4cm x 6.4cm)   | 4       | 10     | 40      |
| 896240 |  | 2.0" x 4.0" (5cm x 10cm)    | 4       | 10     | 40      |
| 896350 |  | 3.0" x 5.0" (7.5cm x 13cm)  | 2       | 10     | 20      |
|        |  | 2.0" x 4.0" (5cm x 10cm)    |         |        |         |
| COM200 |  | 3.0" x 4.0" (7.5cm x 10cm)  | 3       | 10     | 30      |
|        |  | 2.0" x 4.0" (5cm x 10cm)    |         |        |         |

**PALS® Foam** - features a soft foam top.



| PART # | SHAPE                                                                             | DESCRIPTION                | #ELC/PK | #PK/CS | #ELC/CS |
|--------|-----------------------------------------------------------------------------------|----------------------------|---------|--------|---------|
| 975200 |  | 2.0" (5cm)                 | 4       | 10     | 40      |
| 975220 |  | 2.0" x 2.0" (5cm x 5cm)    | 4       | 10     | 40      |
| 975240 |  | 2.0" x 3.5" (5cm x 9cm)    | 4       | 10     | 40      |
| 975250 |  | 2.0" x 5.0" (5cm x 13cm)   | 2       | 10     | 20      |
| 976230 |  | 1.5" x 2.5" (4cm x 6.4cm)  | 4       | 10     | 40      |
| 976240 |  | 2.0" x 4.0" (5cm x 10cm)   | 4       | 10     | 40      |
| 976350 |  | 3.0" x 5.0" (7.5cm x 13cm) | 2       | 10     | 20      |

**PALS® StimTrode®** - features a white polycoated fabric top.

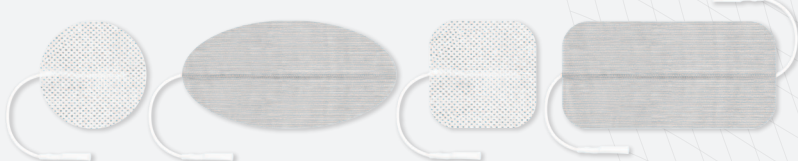

| PART # | SHAPE                                                                               | DESCRIPTION              | #ELC/PK | #PK/CS | #ELC/CS |
|--------|-------------------------------------------------------------------------------------|--------------------------|---------|--------|---------|
| ST32D  |    | 1.25" (3.2cm)            | 4       | 10     | 40      |
| ST50D  |    | 2.0" (5cm)               | 4       | 10     | 40      |
| ST75D  |   | 2.75" (7cm)              | 4       | 10     | 40      |
| ST5050 |  | 2.0" x 2.0" (5cm x 5cm)  | 4       | 10     | 40      |
| ST5090 |  | 2.0" x 3.5" (5cm x 9cm)  | 4       | 10     | 40      |
| ST5010 |  | 2.0" x 4.0" (5cm x 10cm) | 4       | 10     | 40      |
| ST510D |  | 2.0" x 4.0" (5cm x 10cm) | 2       | 10     | 20      |

**PALS® Blue** - specialized gel gentle enough for sensitive skin.




| PART # | SHAPE                                                                               | DESCRIPTION             | #ELC/PK | #PK/CS | #ELC/CS |
|--------|-------------------------------------------------------------------------------------|-------------------------|---------|--------|---------|
| 901220 |  | 2.0" x 2.0" (5cm x 5cm) | 4       | 10     | 40      |
| 901240 |  | 1.5" x 3.5" (4cm x 9cm) | 4       | 10     | 40      |

# UltraStim<sup>®</sup> Electrodes

Patented Current Controlling Neurostimulation

## Experience Stimulation Ease and Comfort

These electrodes are developed with a scientifically calculated silver grid pattern printed on conductive material. This grid pattern precisely controls and optimizes the electrical current over the surface of the electrode. And the border around the pattern allows for an intentional current drop off. This eliminates edge biting and stinging and results in stimulation comfort.

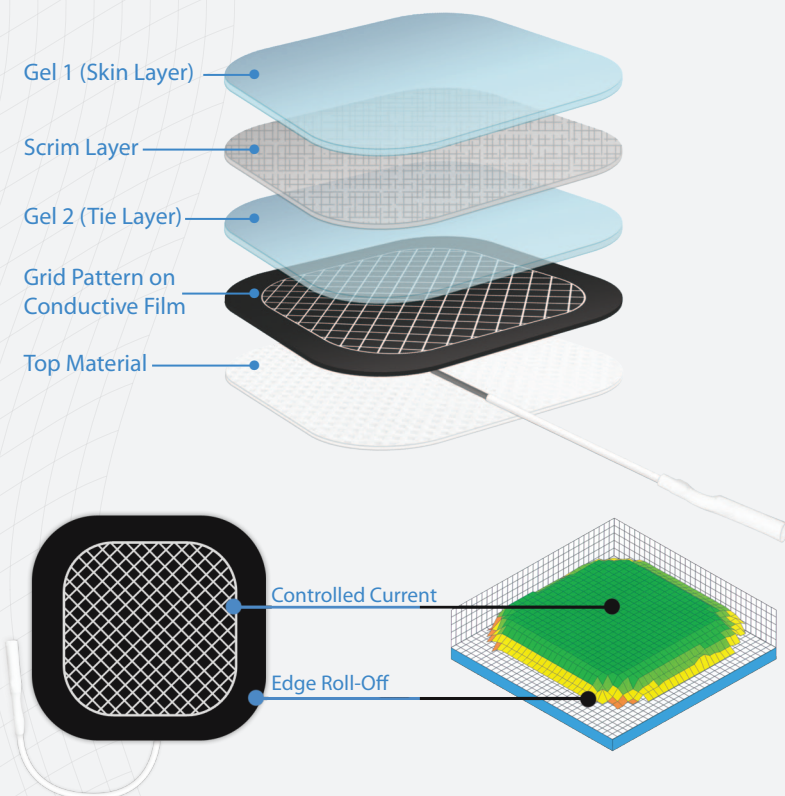

## Featuring MultiStick<sup>®</sup> Hydrogel

Our patented MultiStick hydrogel used in the construction of UltraStim electrodes is specially developed for adherence to conductive film as well as easy release from the patient's skin.

## Available in multiple configurations

UltraStim® PAD




UltraStim® WIRE




UltraStim® SNAP



UltraStim® Universal Dual

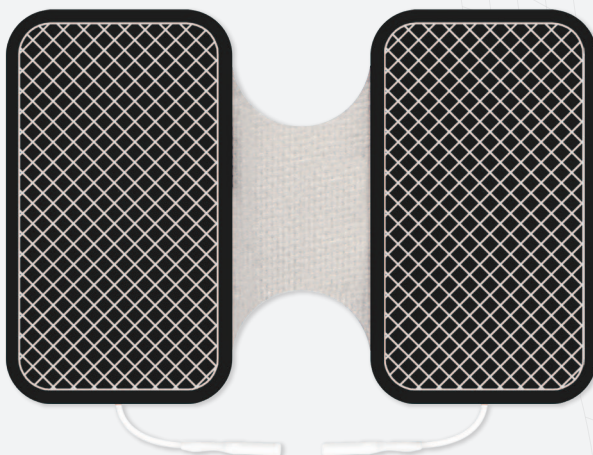

UltraStim® X

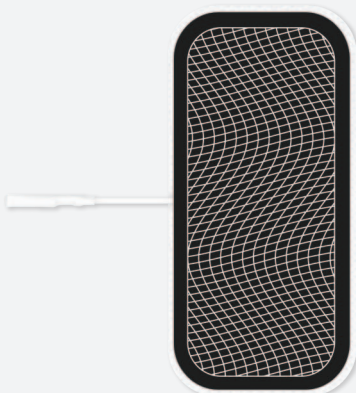


UltraStim® X Blue



# UltraStim® Wire


- Proprietary silver mixture and grid pattern precisely control and optimize electrical current over the electrode surface.
- Areas within the lines of the silver grid pattern are calculated to provide optimum current control.
- Polycoated backing provides an additional moisture barrier.

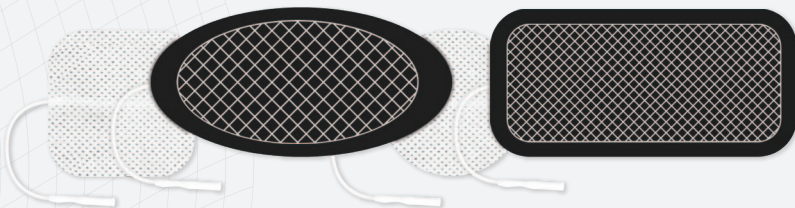

| PART # | SHAPE                                                                             | DESCRIPTION               | #ELC/PK | #PK/CS | #ELC/CS |
|--------|-----------------------------------------------------------------------------------|---------------------------|---------|--------|---------|
| UF1250 |  | 1.25" (3.2cm)             | 4       | 10     | 40      |
| UF2000 |  | 2.0" (5cm)                | 4       | 10     | 40      |
| UF2750 |  | 2.75" (7cm)               | 4       | 10     | 40      |
| UF2020 |  | 2.0" x 2.0" (5cm x 5cm)   | 4       | 10     | 40      |
| UF2040 |  | 2.0" x 4.0" (5cm x 10cm)  | 4       | 10     | 40      |
| UF1525 |  | 1.5" x 2.5" (4cm x 6.4cm) | 4       | 10     | 40      |
| UF2400 |  | 2.0" x 4.0" (5cm x 10cm)  | 4       | 10     | 40      |

# UltraStim® Dual

- Offers a fixed distance between electrodes while distributing current evenly and comfortably.

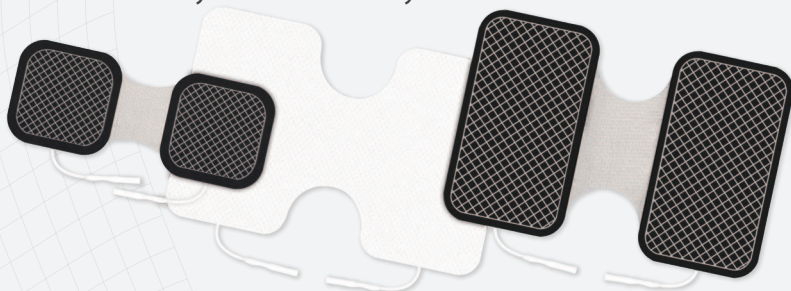


| PART # | SHAPE                                                                               | DESCRIPTION              | #ELC/PK | #PK/CS | #ELC/CS |
|--------|-------------------------------------------------------------------------------------|--------------------------|---------|--------|---------|
| ULB355 |  | 3.5" x 5.0" (9cm x 13cm) | 1       | 10     | 10      |
| ULB250 |  | 2.0" x 5.0" (5cm x 13cm) | 2       | 10     | 20      |

# UltraStim®X

## Patented Neurostimulation Electrodes

### Performance and Value

This electrode features our precisely calculated and proprietary silver grid pattern to control electrical current evenly over the surface of the electrode. The featured oversize top is specially designed to facilitate ease of proper removal to further promote electrode durability.

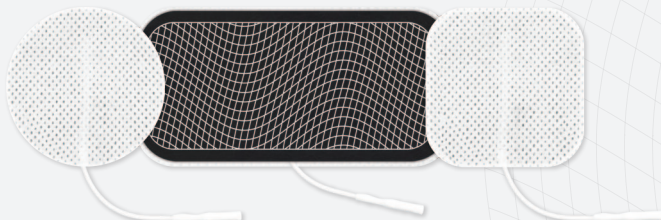


| PART #    | SHAPE                                                                             | DESCRIPTION              | #ELC/PK | #PK/CS | #ELC/CS |
|-----------|-----------------------------------------------------------------------------------|--------------------------|---------|--------|---------|
| USX2000   |  | 2.0" (5cm)               | 4       | 10     | 40      |
| USX2020   |  | 2.0" x 2.0" (5cm x 5cm)  | 4       | 10     | 40      |
| USX2040   |  | 2.0" x 4.0" (5cm x 10cm) | 2       | 10     | 20      |
| USX2040-4 |  | 2.0" x 4.0" (5cm x 10cm) | 4       | 10     | 40      |

# UltraStim®X Blue

## Patented Neurostimulation Electrodes

### Innovation for Sensitive Skin

These electrodes feature BLUE hydrogel for those with sensitive skin. A unique patented silver grid pattern is strategically printed on the top of the electrode, providing a comfortable therapy treatment for patients.

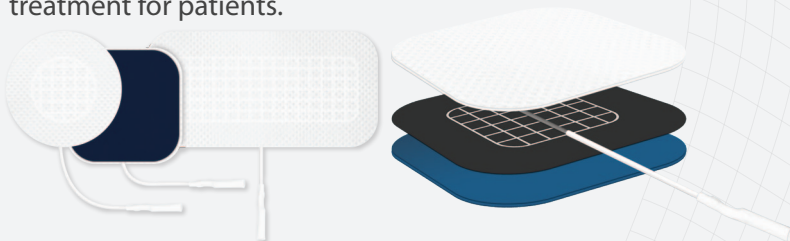

| PART #    | SHAPE                                                                               | DESCRIPTION              | #ELC/PK | #PK/CS | #ELC/CS |
|-----------|-------------------------------------------------------------------------------------|--------------------------|---------|--------|---------|
| USXPB5000 |  | 2.0" (5cm)               | 4       | 10     | 40      |
| USXPB5050 |  | 2.0" x 2.0" (5cm x 5cm)  | 4       | 10     | 40      |
| USXPB5010 |  | 2.0" x 4.0" (5cm x 10cm) | 4       | 10     | 40      |

# UltraStim® Snap

Patented Neurostimulation Electrodes

- Conductive eyelet is scientifically chosen to match impedance in the entire system to provide even current distribution.
- Proprietary silver mixture and grid pattern precisely control and optimize current over the electrode surface.
- Border surrounding the grid pattern allows for an intentional current roll off that eliminates edge biting and stinging.



| PART #   | SHAPE | DESCRIPTION              | #ELC/PK | #PK/CS | #ELC/CS |
|----------|-------|--------------------------|---------|--------|---------|
| SN2020   |       | 2.0" x 2.0" (5cm x 5cm)  | 4       | 10     | 40      |
| SN2040   |       | 2.0" x 4.0" (5cm x 10cm) | 4       | 10     | 40      |
| SN2040-D |       | 2.0" x 4.0" (5cm x 10cm) | 2       | 10     | 20      |


# UltraStim®

Garments & Pad Electrodes

## Easy-to-Use, Patented, Lower Back TENS Solution

- Pad electrodes control current with a scientifically calculated silver grid pattern printed on conductive film and laminated between four layers of hydrogel.
- UltraStim garments are developed specifically for convenient one-person home application. Instructions printed directly on the garment assist in proper and consistent electrode positioning. The electrodes remain on the garment for easy storage and reapplication.

## Four Gel Layers

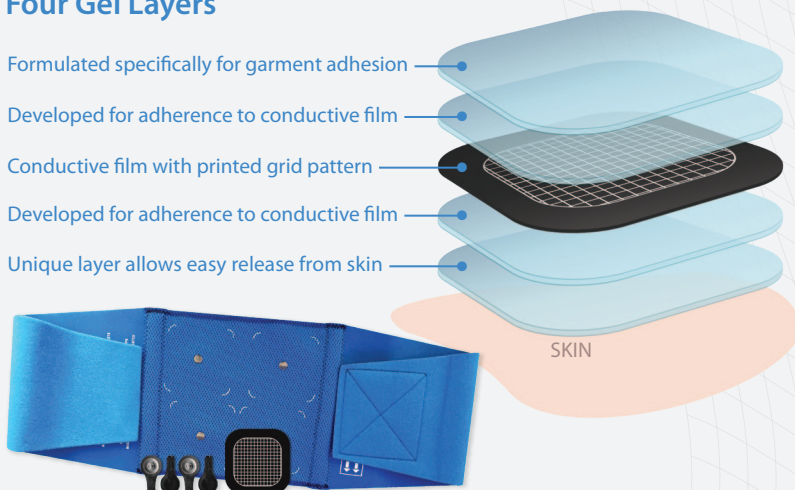

| PART # | DESCRIPTION                                                                                                                                   | PER PK | #PK/CS | #ELC/CS |
|--------|-----------------------------------------------------------------------------------------------------------------------------------------------|--------|--------|---------|
| B22-S  | Back Garment Kit contains:<br>1 S/M Back Garment (Waist 25"-42", 63cm-107cm)<br>4 UltraStim Electrodes 2.0" x 2.0"<br>4 Black Snap Adapters   | 1      | 1      | 1       |
| B22-L  | Back Garment Kit contains:<br>1 L/XL Back Garment (Waist 40"-54", 102cm-137cm)<br>4 UltraStim Electrodes 2.0" x 2.0"<br>4 Black Snap Adapters | 1      | 1      | 1       |
| US2020 | 2.0" x 2.0" (5cm x 5cm)                                                                                                                       | 4      | 10     | 40      |
| SA1004 | Black Snap Adapter                                                                                                                            | 4      | 10     | 40      |
| SA2004 | White Snap Adapter                                                                                                                            | 4      | 10     | 40      |

# ValuTrobe®

## Conductive Film Neurostimulation Electrodes

### Superior Quality at an Economical Price

- Every lead wire interface is customized for each electrode to optimize current distribution in various shapes and sizes.
- Known for their durability and multiple applications to the skin.



#### ValuTrobe® Cloth size and shape options

| PART # | SHAPE | DESCRIPTION                | #ELC/PK | #PK/CS | #ELC/CS |
|--------|-------|----------------------------|---------|--------|---------|
| CF3200 | ○—    | 1.25" (3.2cm)              | 4       | 10     | 40      |
| CF5000 | ○—    | 2.0" (5cm)                 | 4       | 10     | 40      |
| CF7000 | ○—    | 2.75" (7cm)                | 4       | 10     | 40      |
| CF5050 | □—    | 2.0" x 2.0" (5cm x 5cm)    | 4       | 10     | 40      |
| CF4090 | ▭—    | 1.5" x 3.5" (4cm x 9cm)    | 4       | 10     | 40      |
| CF5090 | ▭—    | 2.0" x 3.5" (5cm x 9cm)    | 4       | 10     | 40      |
| CF7515 | ▭—    | 3.0" x 5.0" (7.5cm x 13cm) | 2       | 10     | 20      |
| CF4065 | ▭—    | 1.5" x 2.5" (4cm x 6.4cm)  | 4       | 10     | 40      |
| CF5010 | ▭—    | 2.0" x 4.0" (5cm x 10cm)   | 4       | 10     | 40      |
| CF1020 | ▭—    | 4.0" x 7.0" (10cm x 18cm)  | 1       | 10     | 10      |

#### ValuTrobe® Foam size and shape options

| PART # | SHAPE | DESCRIPTION                | #ELC/PK | #PK/CS | #ELC/CS |
|--------|-------|----------------------------|---------|--------|---------|
| CFF125 | ○—    | 1.25" (3.2cm)              | 4       | 10     | 40      |
| CFF200 | ○—    | 2.0" (5cm)                 | 4       | 10     | 40      |
| CFF275 | ○—    | 2.75" (7cm)                | 4       | 10     | 40      |
| CFF202 | □—    | 2.0" x 2.0" (5cm x 5cm)    | 4       | 10     | 40      |
| CFF153 | ▭—    | 1.5" x 3.5" (4cm x 9cm)    | 4       | 10     | 40      |
| CFF203 | ▭—    | 2.0" x 3.5" (5cm x 9cm)    | 4       | 10     | 40      |
| CFF305 | ▭—    | 3.0" x 5.0" (7.5cm x 13cm) | 2       | 10     | 20      |
| CFF152 | ▭—    | 1.5" x 2.5" (4cm x 6.4cm)  | 4       | 10     | 40      |
| CFF204 | ▭—    | 2.0" x 4.0" (5cm x 10cm)   | 4       | 10     | 40      |

# ValuTrode®*χ*

- These electrodes feature an oversized border to facilitate proper and easy removal, further extending the life of the electrode.



## ValuTrode®*χ* Cloth size and shape options

| PART #   | SHAPE                                                                             | DESCRIPTION              | #ELC/PK | #PK/CS | #ELC/CS |
|----------|-----------------------------------------------------------------------------------|--------------------------|---------|--------|---------|
| VTX5000  |  | 2.0" (5cm)               | 4       | 10     | 40      |
| VTX5050  |  | 2.0" x 2.0" (5cm x 5cm)  | 4       | 10     | 40      |
| VTX50100 |  | 2.0" x 4.0" (5cm x 10cm) | 4       | 10     | 40      |

## ValuTrode®*χ* Foam size and shape options

| PART #    | SHAPE                                                                             | DESCRIPTION              | #ELC/PK | #PK/CS | #ELC/CS |
|-----------|-----------------------------------------------------------------------------------|--------------------------|---------|--------|---------|
| VTX2000   |  | 2.0" (5cm)               | 4       | 10     | 40      |
| VTX2020   |  | 2.0" x 2.0" (5cm x 5cm)  | 4       | 10     | 40      |
| VTXF50100 |  | 2.0" x 4.0" (5cm x 10cm) | 4       | 10     | 40      |

# ValuTrode®*LITE*

- The low profile design maximizes flexibility and allows the electrode to conform to all body contours.
- Feature an anti-lift adhesive safety border that gently adheres to the skin to reduce stinging, edge biting, and roll ups.




| PART # | SHAPE                                                                               | DESCRIPTION              | #ELC/PK | #PK/CS | #ELC/CS |
|--------|-------------------------------------------------------------------------------------|--------------------------|---------|--------|---------|
| VL4500 |  | 2.0" (5cm)               | 4       | 10     | 40      |
| VL4545 |  | 2.0" x 2.0" (5cm x 5cm)  | 4       | 10     | 40      |
| VL4595 |  | 2.0" x 4.0" (5cm x 10cm) | 4       | 10     | 40      |

# EconoStick®

## Patented Adhesive Hydrogel

The thinner version of Axelgaard's world renowned patented MultiStick® hydrogel is the first of its kind. EconoStick®'s advanced, scrimless technology combines two hydrogel layers and results in a thinner hydrogel that is "stickier" to the touch providing optimum adhesion.

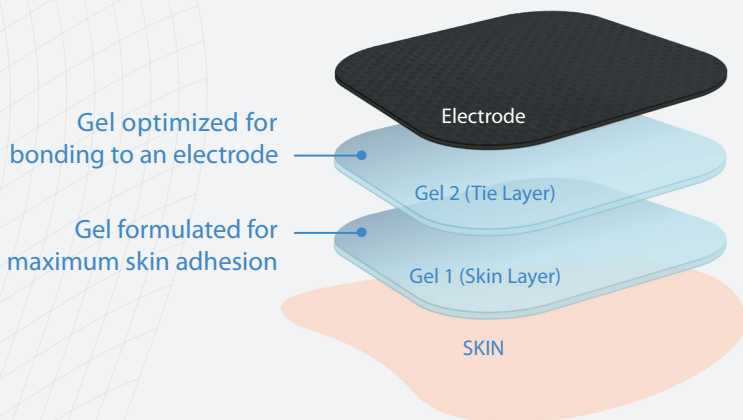

### Scrimless Technology

In addition to optimum adhesion, this scrimless approach results in an electrode with increased flexibility.



# EconoStim®

Low-Profile and Economical

EconoStim® electrodes feature EconoStick® hydrogel. These electrodes include our oversize border to provide easy removal. The EconoStim electrode, with its black top, is the latest addition to our economy electrode product line.



| PART #   | SHAPE                                                                             | DESCRIPTION              | #ELC/PK | #PK/CS | #ELC/CS |
|----------|-----------------------------------------------------------------------------------|--------------------------|---------|--------|---------|
| ESB5000  |  | 2.0" (5cm)               | 8       | 10     | 80      |
| ESB5050  |  | 2.0" x 2.0" (5cm x 5cm)  | 8       | 10     | 80      |
| ESB50100 |  | 2.0" x 4.0" (5cm x 10cm) | 4       | 10     | 40      |

# EconoStim® SELECT

EconoStim® Select with its distinctive black leadwire is well suited for TENS and stimulation of small muscle groups. It is our most economical stimulation electrode.



| PART #  | SHAPE                                                                               | DESCRIPTION             | #ELC/PK | #PK/CS | #ELC/CS |
|---------|-------------------------------------------------------------------------------------|-------------------------|---------|--------|---------|
| ESB3050 |  | 1.2" x 2.0" (3cm x 5cm) | 12      | 10     | 120     |

# **Ionto480<sup>®</sup>** Buffered Iontophoresis | Electrodes



## Iontophoresis Drug Delivery System

- Easily accessible reservoir pad fill area allows for easy application of solution prior to removal of adhesive backing.
- Flexible foam backing with medical grade adhesive conforms to body contours.

## Buffered Return Electrode

- Patented buffered gel minimizes skin irritation.

Kits Offer Active Drug Delivery Electrodes in a Variety of Sizes

| PART # | SHAPE                                                                               | KIT CONTENTS                                                                                                                          | #KITS/CS |
|--------|-------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------|----------|
| DDS275 |  | <b>1 Small Square Active Drug Delivery Electrode</b><br>(Approximate Fill Volume - 1.5cc)<br>1 Buffer Electrode, 2 Alcohol Prep Pads  | 10       |
| DDM325 |  | <b>1 Medium Square Active Drug Delivery Electrode</b><br>(Approximate Fill Volume - 2.5cc)<br>1 Buffer Electrode, 2 Alcohol Prep Pads | 10       |
| DDL350 |  | <b>1 Large Square Active Drug Delivery Electrode</b><br>(Approximate Fill Volume - 4.0cc)<br>1 Buffer Electrode, 2 Alcohol Prep Pads  | 10       |
| ddb375 |  | <b>1 Butterfly Active Drug Delivery Electrode</b><br>(Approximate Fill Volume - 2.0cc)<br>1 Buffer Electrode, 2 Alcohol Prep Pads     | 10       |

1 Instruction, 10 Individually Packaged Kits

# LittlePALS<sup>®</sup>

Our patented neonatal pediatric electrodes are a revolutionary advancement in neonatal apnea, electrocardiogram (ECG), and electromyography (EMG) biofeedback monitoring.

To meet the needs of infant patients, these electrodes provide adhesion, pliability, and sensitivity while offering resistance to heat, humidity, and dry-out conditions.



| PART #  | DESCRIPTION - STANDARD                                                                                          | #ELC/PK | #PK/CS | #ELC/CS |
|---------|-----------------------------------------------------------------------------------------------------------------|---------|--------|---------|
| SEN2001 | <b>0.9" x 1.1" (2.2cm x 2.8cm)</b><br>Green, White, Black 24" Radiotranslucent Wires<br>Socket .059", 1.5mm DIN | 3       | 10     | 30      |
| SEN2002 | <b>0.9" x 1.1" (2.2cm x 2.8cm)</b><br>Red, Yellow, Green 24" Radiotranslucent Wires<br>Socket .059", 1.5mm DIN  | 3       | 10     | 30      |
| SEN4001 | <b>0.9" x 1.1" (2.2cm x 2.8cm)</b><br>Green, White, Black 24" Copper Wires<br>Socket .059", 1.5mm DIN           | 3       | 10     | 30      |
| SEN5001 | <b>0.9" x 1.1" (2.2cm x 2.8cm)</b><br>Snap, White Fabric Top                                                    | 6       | 10     | 60      |

| PART #  | DESCRIPTION - DUAL SENSOR                                                                   | #ELC/PK | #PK/CS | #ELC/CS |
|---------|---------------------------------------------------------------------------------------------|---------|--------|---------|
| SEN1003 | <b>0.9" x 1.1" (2.2cm x 2.8cm)</b><br>2 black, 2 white 24" Copper Wire<br>Socket .080", 2mm | 4       | 10     | 40      |
| SEN3001 | <b>0.9" x 1.1" (2.2cm x 2.8cm)</b><br>Snap, White Fabric Top                                | 6       | 10     | 60      |



Headquarters: Fallbrook, CA 92028 USA | Europe: DK-8520 Lystrup, Denmark  
Phone: +1 (760) 451-8000 | Phone: +45 8674 2088

Email: [support@axelgaard.com](mailto:support@axelgaard.com)

Website: [www.axelgaard.com](http://www.axelgaard.com)



Intertek

Axelgaard Logo, MultiStick, PALS, UltraStim, ValuTrove, StimTrove, EconoStim, EconoStick and Ionto480 are Registered Trademarks of Axelgaard Manufacturing Co., Ltd.

